

LAB

CRP Ultra


One product for ultra sensitive
and standard procedure



SENTINEL
DIAGNOSTICS

CRP Ultra


CRP Ultra, Inter Assay Precision

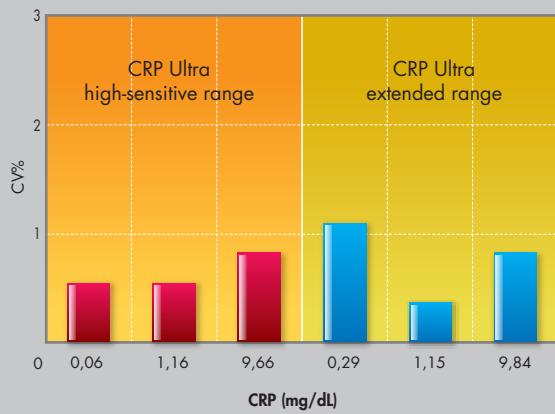

Introduction

C-reactive protein is a non-specific, acute-phase reactant used to diagnose bacterial infectious diseases and inflammatory disorders such as acute rheumatic fever and rheumatoid arthritis. Recent development of a "high sensitivity" assay for CRP has enabled accurate assays at even low levels; the C-reactive protein level is a strong predictor of a coronary event.

CRP Ultra, Risk assesment

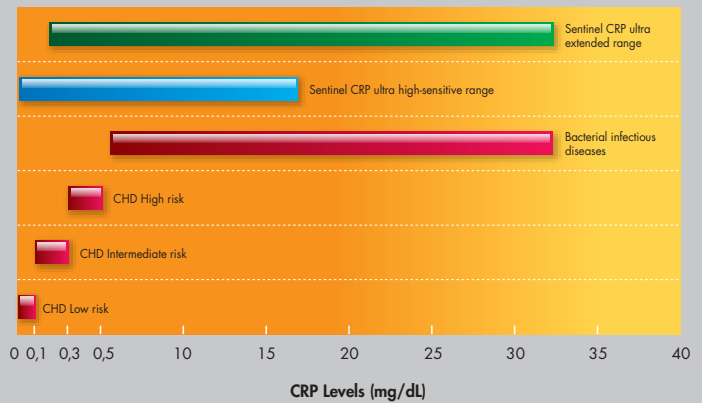

Features and Advantages

- Two different procedure high-sensitive range and extended range to satisfy any user needs by simply changing the calibration points and the sample volume
- Excellent method precision
- No prozone risk up to 100 mg/dL
- No interference up to: RF 550 IU/mL; Bilirubin 30 mg/dL; Hemoglobin 500 mg/dL; Triglycerides 1500 mg/dL
- Prolonged calibration stability, up to 30 days
- Extended reagent on board stability, up to 60 days
- Filled in Universal Vials* suitable for AU series, ADVIA series, Hitachi Modular, Hitachi 917, Dirui series

Method

Sentinel CRP Ultra is a particle-enhanced immunoturbidimetric assay.

Ref	Description	Kit Size
11508H*	CRP Ultra	2x20 mL; 2x20 mL
11508A	CRP Ultra	3x50 mL; 3x50 mL
11532	CRP Calibrator set (0.5-1-2-4-8-16-32 mg/dL)	7x2 mL
11533	CRP Calibrator US (0,25 mg/dL)	1x2 mL
11210	Immuno Control 1	1x5 mL
112121	Immuno Control 2	1x5 mL
11211A	Immuno Control Set (2 Levels)	2x(2x1) mL
11540	CRP Control US (0,05 mg/dL)	3x2 mL

Sentinel CH. SpA

Via Robert Koch, 2
20152 Milan - Italy
Phone: +39 02 3455141
Fax: +39 02 34551464
sentinel@sentinel.it
www.sentinel diagnostics.com



WATCHING OVER LIFE